


Performing by design

Lynx

Rx PTA Balloon Catheter

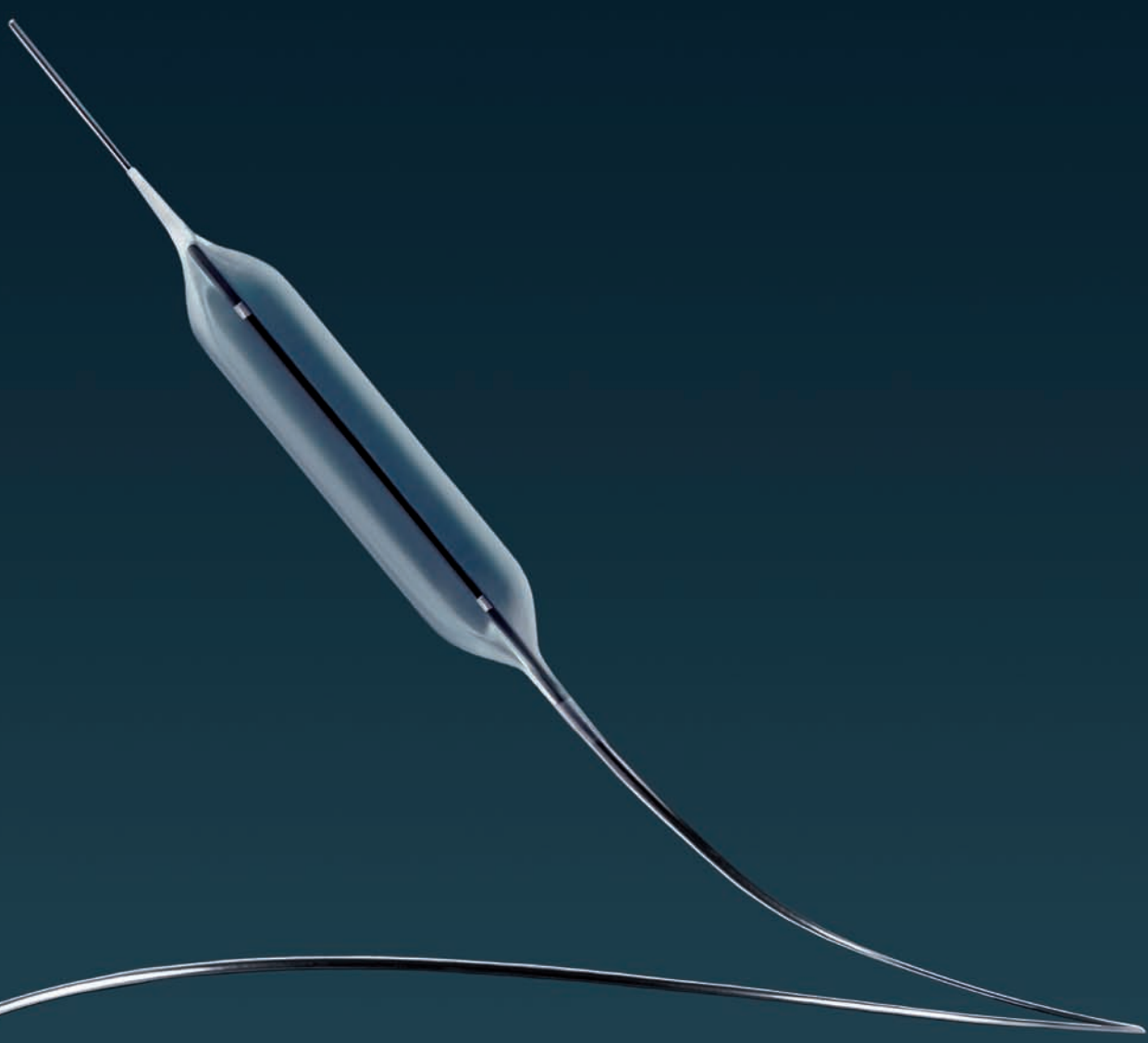

Your best choice
for peripheral procedure

ARTHESYS

Committed to excellence in developing Arterial Therapeutic Systems

Technical specifications

Guidewire compatibility

0.014" (OD ≤ 0.38mm) maximum recommended

Guiding catheter compatibility

- For 5.0 balloon: 0.058" (1.47 mm)
- For 6.0 balloon: 0.071" (1.80 mm)
- For 7.0 balloon: 0.081" (2.06 mm)



BALLOON NOMINAL LENGTH	Balloon nominal diameter		
	5.0 mm	6.0 mm	7.0 mm
10 mm	□	□	□
15 mm	□	□	□
20 mm	□	□	□
25 mm	□	□	□
30 mm	□	□	□
35 mm	□	□	□
40 mm	□	□	□

□ available

PRESSURE (BARS)	Balloon diameter		
	5.0 mm	6.0 mm	7.0 mm
2	4.80	5.76	6.72
3	4.85	5.82	6.79
4	4.90	5.88	6.86
5	4.95	5.94	6.93
6*	5.00	6.00	7.00
7	5.04	6.05	7.05
8	5.08	6.09	7.11
9	5.11	6.14	7.16
10	5.15	6.18	7.21
11	5.19	6.23	7.26
12	5.23	6.27	7.32
13	5.26	6.32	7.37
14**	5.30	6.36	7.42
15	5.34	6.41	7.47
16	5.38	6.45	7.53

* Nominal Pressure

** Rated Burst Pressure

Optimal trackability

Innovative hydrophilic coating

Excellent pushability

Hypotube designed with specific treatment

Hypotube shaft

Especially devised to resist the kinking and increase the pushability

6 folding balloon

For an easier catheter withdrawal and homogeneous stent deployment

